

KIT COMPONENTS

Part Number	Description	Qty.
029-2773-XXX	LH Wrist Rest Assembly	1
029-2776-XXX	RH Wrist Rest Assembly	1
003-1284-00	I / O Instructions.....	1

REQUIRED TOOLS

Thin blade screwdriver
3/8" Wrench


INTRODUCTION

To provide support for wrists during examination and general procedures.

NOTE

9A311001 Wrist Rest Kit can only be used on Models 419-003, 004, 005, and 006 with the new triple articulating headrest.

INSTALLATION

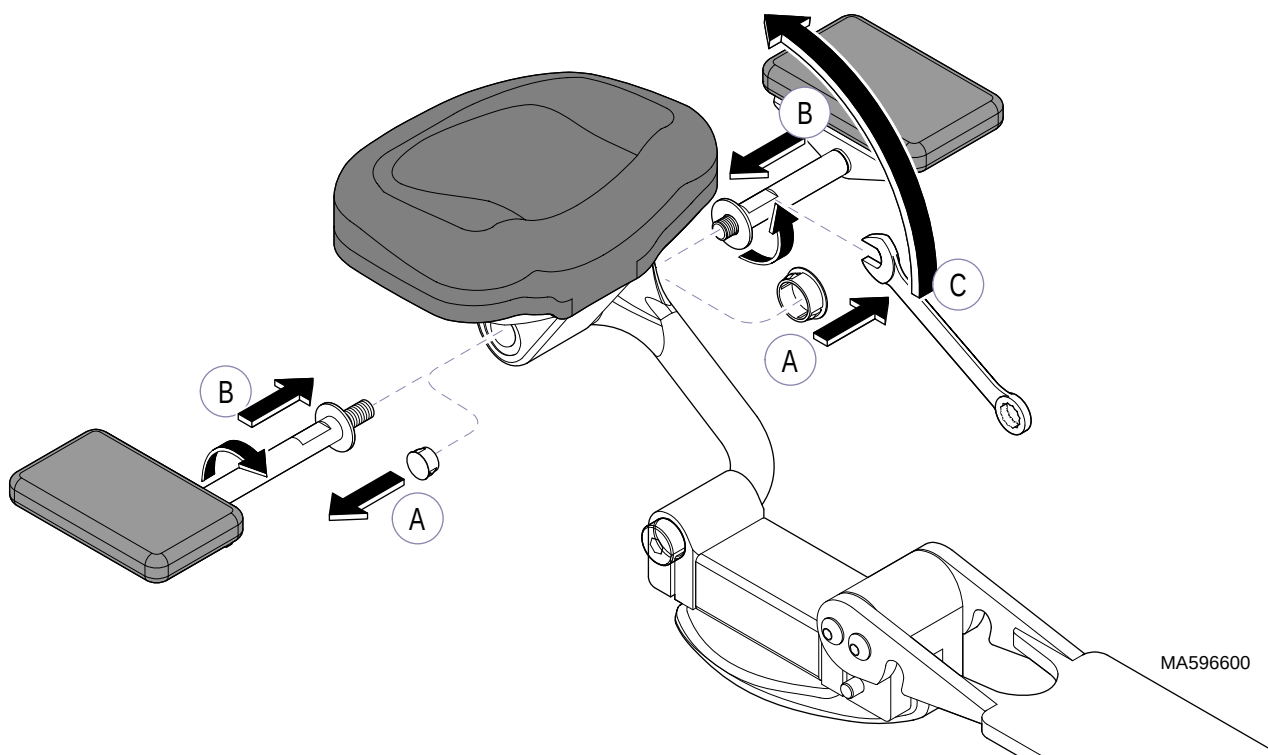

- Remove hole plugs from head rest.
- Hand thread each wrist rest assembly into mounting holes.
- Tighten wrist rest assembly using a 3/8" wrench.